



Laboratory In-Line Mixer

Verso

Capacity

Self-pumping throughputs of up to 7 gallons per minute.

No Bypassing

Product cannot bypass the precision rotor/stator workhead.

Interchangeable Workheads

Single and multistage rotor/stator configurations as standard.

Materials of Construction

Product contact parts in 316L stainless steel.

Motor Specification

1hp (750W) 110 volt, single phase (220 volt optional), 60 Hz. Nominal maximum speed 9000 rpm.

Speed Control


Infinitely variable electronic speed control with integral on/off switch.

Inlet and Outlet Connections

1/2" Tri-Clamp connections. Self-draining tangential outlet.

Laboratory In-Line Mixer

Verso

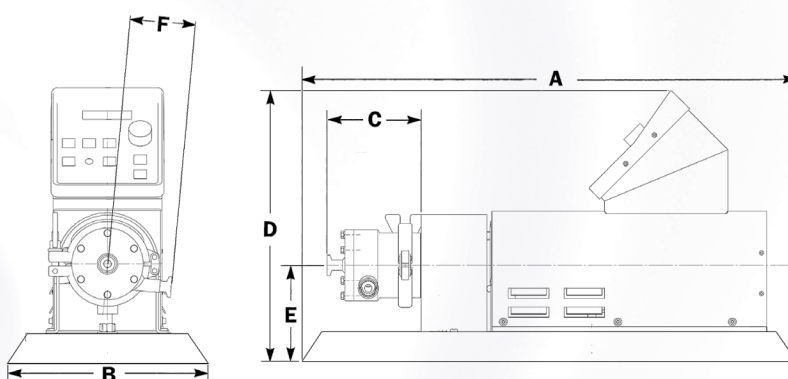

hp	RPM	A	B	C	D	E	F
1	8900	23 ⁵ / ₈	7 ⁷ / ₈	4 ⁵ / ₈	13 ⁵ / ₈	5	3 ³ / ₁₆

All dimensions shown are in inches. The dimensions shown are approximate only and certified diagrams should be used for installation purposes. Silverson reserves the right to change dimensions and specifications without notice.